

# Econo Traceable® Refrigerator Thermometer Instructions

## SPECIFICATIONS

**Display:** 1¼"-High LCD (½" digits)

**Accuracy:** ± 1°C (0 to 25°C); otherwise ± 1.5°C

**Resolution:** 0.1°F or °C

### Temperature

**Range:** 14.2 to 131.0° F (-9.9 to 55.0°C)


**Readings:** display updates every 30 seconds

**Function:** temperature

**Accessories:** Always-Hold™ suction cup, extra battery, Velcro®

**Size:** 1 x 2½ x ½"

**Weight:** ½ oz.

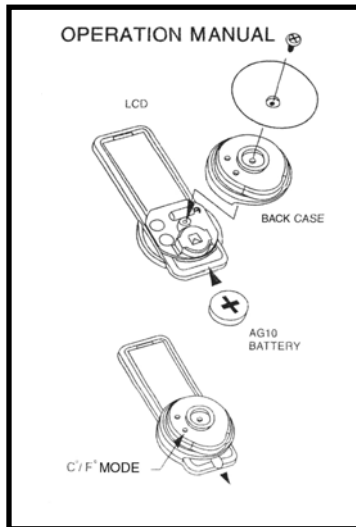

### Note: Before Using

1. Unit is shipped with battery installed. Unit is ON and displaying current ambient temperature.
2. If display is not showing temperature, replace battery with extra one supplied. (See Battery Replacement.)

### Change Reading to Fahrenheit (F°) or Celsius (C°): MODE

1. Hold unit in one hand with back facing up. Slightly fold back suction cup. MODE opening is visible on lower left side of unit back.
2. With other hand, insert tip of bent paper clip in MODE opening. Press paper clip firmly and hold for one second. Remove paper clip. Reading will change from F° to C° or C° to F°.

**Note:** SET opening on back of unit has no function on this model.

### Mounting Unit

**Note:** Because refrigerators have varying internal air temperatures, unit should be placed as close as possible to materials being observed.

1. Free-standing – unit may be set on any flat surface.
2. Suction cup – wet suction cup with warm water. Apply to any clean surface.
3. Velcro® – remove paper cover from one piece of Velcro®. Press Velcro® firmly onto suction cup. (Suction cup will flatten.) Remove tape from other Velcro® piece. Press Velcro® firmly onto selected surface for mounting unit. Press unit firmly in place.

**Note:** Temperatures above 100°F are indicated by “100” at the top of the display.

### All Operational Difficulties

If this unit does not function properly for any reason, replace battery with fresh 1.5 volt silver-oxide battery. (See Battery Replacement.) Low battery power can occasionally cause any number of “apparent” operational difficulties. Replacing battery with a fresh one solves most difficulties. Equivalent battery replacements are 1.5 volt silver-oxide battery. Replacement battery Cat. No. 1105.

**Note:** If unit is exposed to rapid changes in temperature, moisture may condense on internal electronic parts or on display. (Symptoms include erratic, incorrect, or no readings on display.) To dry out, open unit by following steps 1,2 in Battery Replacement. Allow to dry at room temperature for 24 hours. Close unit by following steps 6,7 in Battery Replacement.

**Note:** If unit is operating in temperatures below recommended range, liquid crystal display may freeze.

### Battery Replacement

1. To open unit, grasp suction cup and turn counter-clockwise until it can be removed or use a screwdriver to remove the screw.
2. Remove BACK CASE.
3. Hold unit in one hand.
4. With other hand, use tip of paper clip to put pressure on top of battery. Press firmly on top of battery until it slides out.
5. With +side facing up, slide fresh battery in place.
6. Replace BACK CASE.
7. Grasp suction cup. Replace screw in opening in center of BACK CASE. Turn suction cup by hand, clockwise or use a screwdriver until it is firmly in place.

**Cat. No. 4157**

**Econo Traceable® Refrigerator Thermometer Instructions**

**92-4157-00**

Rev. 0 083005

For more information, please contact us:

[ExpotechUSA](#)  
[10700 Rockley Road](#)  
[Houston, Texas 77099](#)  
[USA](#)

[281-496-0900 \[voice\]](#)

[281-496-0400 \[fax\]](#)

E-mail: [sales@expotechusa.com](mailto:sales@expotechusa.com)

Website: [www.ExpotechUSA.com](http://www.ExpotechUSA.com)