

For more information, please contact us:

[ExpotechUSA](#)
[10700 Rockley Road](#)
[Houston, Texas 77099](#)
[USA](#)

[281-496-0900 \[voice\]](#)

[281-496-0400 \[fax\]](#)


E-mail: sales@expotechusa.com

Website: www.ExpotechUSA.com

The Orion basic portable pH meters offer fast and accurate measurements at a value price. All meter functions are identified by color coded keys allowing easy operation.

Orion 210Aplus

pH/Temperature

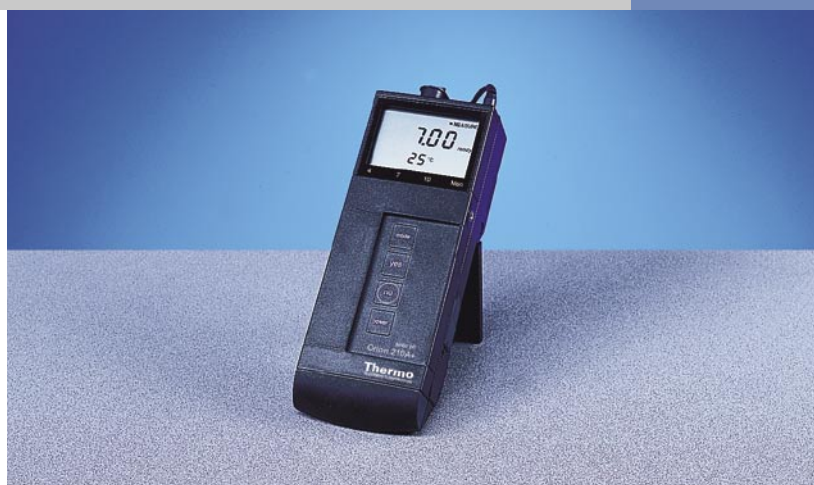

Basic Portable pH/Temperature Meter Orion 210Aplus

- 2 point pH Calibration
- Simple operation yet accurate results
- pH Resolution to 0.01

The Orion 210Aplus is an affordable, durable portable meter for performing routine pH/temperature measurements in the field, plant or laboratory. Housed in a splash- and chemical-resistant case, this meter includes a built-in electrode storage sleeve and cable management. Simple four-button operation makes this the meter of choice for many applications. The two step calibration with built-in

Auto-Buffer-Recognition automatically calculates and displays electrode slope value. This gives the user valuable electrode status feedback to ensure quality results. Meter and electrode package includes a Low Maintenance™ pH electrode. Three year limited warranty.

Orion 210Aplus

Orion 210Aplus Specifications

		Specification	Benefit
pH	Range	0.00 - 14.00	pH in any sample
	Resolution	0.01/0.1	
	Relative Accuracy	± 0.02	
	Slope	80 to 120%	Ensure quality of results
	Auto-Buffer-Recognition	4.01, 7.00, 10.01	Quicker calibrations bracket your sample's pH
	Calibration Points	2	Improved accuracy
Temperature	Range	-5 to 105 °C	
	Resolution	1 °C	
	Relative Accuracy	1 °C	
E_H ORP mode		No	
Display		Custom LCD	Easy to read
Inputs		One BNC, one pin tip ref, ATC, power	
Power Requirements		One 9 V alkaline/lithium battery, or line adapter for 110 V, 220 V, or 240 V	
Agency Approvals		CE, CSA, UL	World wide operation
Environmental Requirements		5 to 45 °C and 5 to 85% relative humidity, non-condensing	

Above data generated under optimal/controlled circumstances.

Orion 210Aplus Ordering Information

Package Includes	Cat. No. for Complete Meter Package	Electrode Cat. No.
Meter and Electrode	0210A0 with epoxy body pH electrode*	9106BNWP
Meter Only	0210A2*	

* Replaceable battery, 9 V included.

All Meters and Meter Kits must be purchased with a Plug/Wall adapter at the proper voltage

Adapter		Cat. No.
Plug/Wall adapter, 110 V	Order separately	020125
Plug/Wall adapter, 220 V	Order separately	020130
Plug/Wall adapter, 240 V	Order separately	020135